

# Subretinal Photovoltaic Implant PRIMA Provides Central Vision in Subjects with Geographic Atrophy Due to Age-related Macular Degeneration

José-Alain Sahel,<sup>1,2,3</sup> Yannick Le Mer,<sup>2,3</sup> Frank G. Holz,<sup>4</sup> Joseph Martel,<sup>1</sup> Lisa Olmos de Koo,<sup>5</sup> Mahi Muqit<sup>6,7</sup> and Daniel Palanker<sup>8</sup>

1. Department of Ophthalmology, University of Pittsburgh School of Medicine, Pittsburgh, PA, USA; 2. Department of Ophthalmology, Fondation Ophthalmologique A. de Rothschild, Paris, France; 3. Clinical Investigation Center, Quinze-Vingts National Eye Hospital, Paris, France; 4. Department of Ophthalmology, University of Bonn, Bonn, Germany; 5. Department of Ophthalmology, University of Washington, Seattle, WA, USA; 6. Vitreoretinal Service, Moorfields Eye Hospital, London, UK; 7. Institute of Ophthalmology, University College London, London, UK; 8. Department of Ophthalmology, Stanford University, Stanford, CA, USA

**A**ge-related macular degeneration is one of the leading causes of irreversible vision loss in ageing populations worldwide. The PRIMA bionic vision system is a photovoltaic substitute for lost photoreceptors that is designed for patients with atrophic dry age-related macular degeneration. The subretinal implant, activated by near-infrared light from augmented-reality glasses, provides central prosthetic vision, perceived simultaneously with peripheral natural vision. This editorial reviews PRIMA's mode of action, clinical trial data and potential impact on clinical practice.

## Keywords

Age-related macular degeneration, photoreceptors, photovoltaic implant, retinal prosthesis, subretinal surgery, vision rehabilitation, vision restoration

**Disclosures:** José-Alain Sahel is co-founder of Pixium Vision, GenSight Biologics, Prophesee, Tilak Healthcare, Chronolife, Sparing Vision, Avista Rx, Vegavect, Tenpoint, Cilensee, and SharpEye and Consultant for Tenpoint and Avista Rx. Yannick Le Mer is a Consultant for Pixium Vision. Daniel Palanker is a Consultant for and receives Patent Royalties from Pixium Vision. Frank Holz, Joseph Martel, Lisa C Olmos de Koo and Mahi Muqit have no financial or non-financial relationships or activities to declare in relation to this article.

**Review process:** Double-blind peer review.

**Compliance with ethics:** This article is an opinion piece and does not report on new clinical data, or any studies with human or animal subjects performed by any of the authors.

**Data availability:** Data sharing is not applicable to this article as no datasets were generated or analysed during the writing of this article.

**Authorship:** The named authors meet the International Committee of Medical Journal Editors (ICMJE) criteria for authorship of this manuscript, take responsibility for the integrity of the work as a whole, and have given final approval for the version to be published.

**Access:** This article is freely accessible at touchOPHTHALMOLOGY.com. © Touch Medical Media 2024.

**Received:** 7 July 2023

**Accepted:** 25 September 2023

**Published online:** 30 January 2024

**Citation:** touchREVIEWS in Ophthalmology. 2024;18(1):Online ahead of journal publication

**Corresponding author:** José-Alain Sahel, Department of Ophthalmology, University of Pittsburgh School of Medicine, 1622 Locust St., Pittsburgh, PA, 15219, USA. E: j.sahel@gmail.com


**Support:** No funding was received in the publication of this article.

Age-related macular degeneration (AMD) is one of the main causes of irreversible vision loss in ageing populations worldwide.<sup>1</sup> In 2019, the US recorded 1.49 million people aged 40 years or older with late-stage AMD, reflecting a crude prevalence rate of 0.94%.<sup>2</sup> Geographic atrophy (GA), the late stage of non-neovascular AMD, is characterized by gradual and progressive atrophy of photoreceptors, while most of the inner retinal cells remain intact.<sup>3</sup> The central atrophic area expands over time, leading to increased vision loss and a reduced ability to read and recognize faces, severely impacting patients' quality of life.<sup>4</sup> The PRIMA bionic vision system (Pixium Vision SA, Paris, France) is a photovoltaic wireless subretinal implant designed to improve central vision in patients with atrophic AMD.<sup>5</sup> Currently, PRIMA is the only experimental therapy capable of restoring vision in patients with degenerated photoreceptors, such as those affected by atrophic AMD. This editorial provides an overview of (1) current treatment strategies for AMD, (2) PRIMA's mode of action, clinical trial data and future steps, and (3) its potential impact on clinical practice.

Treatment strategies for AMD should aim at restoring functional central vision without jeopardizing the surrounding retina while also allowing for their simultaneous use.<sup>4</sup> In this context, electronic retinal prostheses enable the reintroduction of visual information into the degenerate retina by electrical stimulation of the remaining neurons. Electrode arrays can be placed either epiretinally to target the retinal ganglion cells or subretinally to stimulate bipolar cells, which comprise the first layer of neurons above the photoreceptors.<sup>5</sup> Argus® II Retinal Prosthesis System (Second Sight Inc, Sylmar, CA, USA), an epiretinal implant approved for patients blinded by retinitis pigmentosa (RP), contains 60 electrodes positioned relatively distant from the target retinal ganglion cells.<sup>6,7</sup> The subretinal implant Alpha IMS and Alpha AMS (Retina Implant AG, Reutlingen, Germany) contain up to 1,600 pixels and has been tested in patients with RP.<sup>8,9</sup> Visual acuity measured using Landolt C-rings was measurable in 4 out of 29 patients with an Alpha IMS, ranging from 20/2000 to 20/546 and in 2 out of 15 patients with an Alpha AMS, ranging from 20/1111 to 20/546.<sup>8,9</sup>

The PRIMA bionic vision system is a wireless subretinal prosthesis in which photovoltaic pixels directly convert projected light patterns into local electric current. It aims to partially restore central vision in subjects affected by GA due to AMD. In addition to the subretinal implant, it includes a pair of glasses equipped with a camera and a digital projector (Figure 1) and a pocket computer for image processing.<sup>4</sup> Images captured by the camera and simplified by the pocket processor are projected from the glasses onto an 18-degree-wide central field on the retina using pulsed near-infrared light (880 nm). The subretinal implant substitutes the photoreceptors lost in GA with photovoltaic pixels, which convert the pulsed near-infrared light into electric current to excite the second-order neurons (bipolar cells) in the retina, providing formed visual perception. The first generation of the PRIMA implant is 2 mm in width and 30 µm in thickness and contains 378 hexagonal pixels, each with a size of 100 µm.<sup>4</sup>

Figure 1: The PRIMA bionic vision system for providing central vision to subjects with geographic atrophy. (a) PRIMA glasses with the projection module including the camera; (b) artistic rendering of the augmented reality glasses capturing a video of a flower and projecting it onto the retina using near-infrared laser.

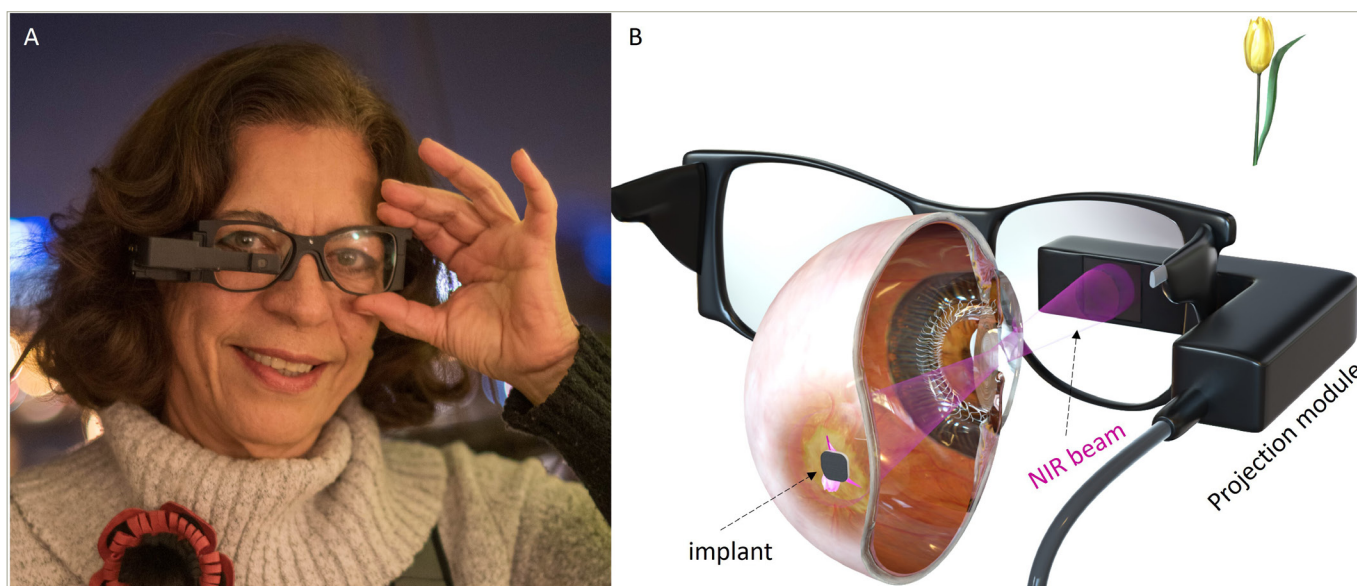

Reproduced and adapted (cropped) with permission from Palanker et al.<sup>4</sup> under CC BY-NC 4.0 DEED (<http://creativecommons.org/licenses/by-nc/4.0>).

A two-phase study evaluating the PRIMA bionic vision system was conducted on five subjects with a GA exceeding three optic disc diameters (ClinicalTrials.gov identifier: NCT03333954). Implantation lasted approximately 2 h, and the subjects' natural visual acuity in the operated eye did not decrease after the surgery.<sup>5</sup> The first phase started after surgery and lasted up to 18 months. Initially, subjects used virtual reality glasses (PRIMA-1), blocking the natural residual vision. The primary endpoint, the perceptual threshold elicited by the implant, was measured in all five subjects.<sup>5</sup> In three subjects with the implant centrally placed in the subretinal space, prosthetic visual acuity closely matched the pixel size, ranging from 20/460 to 20/550 (between 1.1 and 1.3 pixels).<sup>5</sup> The second phase lasted from month 18 to month 36. Subjects used augmented reality glasses (PRIMA-2; Figure 1), allowing the simultaneous use of prosthetic vision and residual natural vision by the fellow eye and the operated eye.<sup>4</sup> Using electronic magnification, the subjects' prosthetic acuity was significantly higher than the residual natural acuity: they were able to recognize optotypes with acuity between 20/98 and 20/63 and to read words that they could not normally read without the prosthesis.<sup>4</sup> The highest improvement in visual acuity was of logMAR 0.78 (20/71 with PRIMA-2 and electronic magnification versus 20/428 without the system) reported at 18–24 months. This feasibility study demonstrated promising results, showcasing the safety, stability and efficacy of the PRIMA subretinal implant over 36 months.<sup>4</sup>

Compared with previous trials using subretinal implants,<sup>9,10</sup> PRIMA has shown better and more consistent performance, which might be partially attributable to patient selection. In patients with AMD, the inner retina in the central scotoma may be better preserved than at the end stage of RP, which is commonly associated with a nearly complete loss of photoreceptors in the entire retina and associated retinal rewiring.<sup>5</sup> Patients using the Argus II system have been reported to have lower visual acuity, measuring at 20/1260 or worse.<sup>11</sup> This implant not only activates proximal neurons but may also stimulate the axons of remote retinal ganglion cells passing underneath the electrodes, causing distorted visual percepts.<sup>12</sup> With PRIMA, prosthetic

visual acuity was found to closely match the pixel size, with an average value of 20/500 and 1.2 pixels, respectively), and the use of electronic magnification enabled optotype recognition, achieving an equivalent acuity superior to 20/100.<sup>4</sup> Furthermore, the epiretinal Argus II and the subretinal Alpha IMS/AMS implants are powered via a transscleral cable, which requires complex surgical procedures and has been associated with postoperative complications.<sup>8,11</sup> In patients with AMD, such subretinal cable may affect the residual photoreceptors and associated peripheral vision. The wireless PRIMA implant simplifies the surgical procedures, which can be performed within only 2 h instead of 4 h with the Argus II system<sup>9,11</sup> and 6–8 h with the Alpha IMS/AMS implants.<sup>8</sup>

Recent studies on pharmacological,<sup>13</sup> gene<sup>14</sup> and stem cell<sup>15</sup> therapies in AMD have focused on slowing or preventing cell loss in retinal degeneration, and there is no established therapy to restore vision in atrophic AMD to date. Compared with other retinal prostheses,<sup>16–18</sup> PRIMA provides better performance and reading capability. Currently, PRIMA is the only therapy to restore vision in a degenerated area of the retina, and it has the potential to become a standard therapy for advanced atrophic AMD. Recently, the US Food and Drug Administration has granted the Breakthrough Device Designation to the PRIMA system.<sup>19</sup> This designation will help Pixium Vision to expedite its development and work closely with the US Food and Drug Administration on refining its regulatory submission in the US. If the system is approved, this designation will also allow better accessibility to certain outpatient and inpatient reimbursement pathways. The PRIMA system is currently under clinical testing in the multicentre PRIMAvra study to demonstrate its performance and confirm its safety in 38 patients with atrophic AMD (ClinicalTrials.gov identifier: NCT04676854).<sup>20</sup> Furthermore, second-generation PRIMA implants, which are now being optimized for future clinical trials, may enhance prosthetic visual acuity up to 20/100 and, with electronic magnification, potentially up to 20/20.<sup>21</sup> Such prosthesis would provide very significant improvement of central vision for a large number of patients blinded by atrophic AMD. □

1. Wong WL, Su X, Li X, et al. Global prevalence of age-related macular degeneration and disease burden projection for 2020 and 2040: A systematic review and meta-analysis. *Lancet Glob Health*. 2014;2:e106–16. DOI: 10.1016/S2214-109X(13)70145-1.
2. Rein DB, Wittenborn JS, Burke-Conte Z, et al. Prevalence of age-related macular degeneration in the US in 2019. *JAMA Ophthalmol*. 2022;140:1202–8. DOI: 10.1001/jamaophthalmol.2022.4401.
3. Fleckenstein M, Mitchell P, Freund KB, et al. The progression of geographic atrophy secondary to age-related macular degeneration. *Ophthalmology*. 2018;125:369–90. DOI: 10.1016/j.ophtha.2017.08.038.
4. Palanker D, Le Mer Y, Mohand-Said S, Sahel JA. Simultaneous perception of prosthetic and natural vision in AMD patients. *Nat Commun*. 2022;13:513. DOI: 10.1038/s41467-022-28125-x.
5. Palanker D, Le Mer Y, Mohand-Said S, et al. Photovoltaic restoration of central vision in atrophic age-related macular degeneration. *Ophthalmology*. 2020;127:1097–104. DOI: 10.1016/j.ophtha.2020.02.024.
6. Luo YH-L, da Cruz L. The Argus® II Retinal Prosthesis System. *Prog Retin Eye Res*. 2016;50:89–107. DOI: 10.1016/j.preteyeres.2015.09.003.
7. Ahuja AK, Yeoh J, Dorn JD, et al. Factors affecting perceptual threshold in Argus II retinal prosthesis subjects. *Transl Vis Sci Technol*. 2013;2:1. DOI: 10.1167/tvst.2.4.1.
8. Stingl K, Bartz-Schmidt KU, Besch D, et al. Subretinal visual implant Alpha IMS – Clinical trial interim report. *Vision Res*. 2015;111:149–60. DOI: 10.1016/j.visres.2015.03.001.
9. Stingl K, Schippert R, Bartz-Schmidt KU, et al. Interim results of a multicenter trial with the new electronic subretinal implant Alpha AMS in 15 patients blind from inherited retinal degenerations. *Front Neurosci*. 2017;11:445. DOI: 10.3389/fnins.2017.00445.
10. Daschner R, Rothermel A, Rudolf R, et al. Functionality and performance of the subretinal implant chip Alpha AMS. *Sensor Mater*. 2018;30:179. DOI: 10.18494/SAM.2018.1726.
11. Humayun MS, Dorn JD, da Cruz L, et al. Interim results from the international trial of second sight's visual prosthesis. *Ophthalmology*. 2012;119:779–88. DOI: 10.1016/j.ophtha.2011.09.028.
12. Nanduri D, Fine I, Horsager A, et al. Frequency and amplitude modulation have different effects on the percepts elicited by retinal stimulation. *Invest Ophthalmol Vis Sci*. 2012;53:205–14. DOI: 10.1167/iovs.11-8401.
13. Shughoury A, Sevgi DD, Ciulla TA. The complement system: A novel therapeutic target for age-related macular degeneration. *Expert Opin Pharmacother*. 2023. DOI: 10.1080/14656566.2023.2257604.
14. American Academy of Ophthalmology. Gene therapy for AMD. 2022. Available at: [www.aao.org/eyenet/article/gene-therapy-for-amd](http://www.aao.org/eyenet/article/gene-therapy-for-amd) (Date last accessed: 23 November 2023).
15. Liu Q, Liu J, Higuchi A. hPSC-derived RPE transplantation for the treatment of macular degeneration. *Prog Mol Biol Transl Sci*. 2023;199:227–69. DOI: 10.1016/bs.pmbts.2023.02.010.
16. Ho AC, Humayun MS, Dorn JD, et al. Long-term results from an epiretinal prosthesis to restore sight to the blind. *Ophthalmology*. 2015;122:1547–54. DOI: 10.1016/j.ophtha.2015.04.032.
17. Daschner R, Rothermel A, Rudolf R, et al. Functionality and performance of the subretinal implant chip alpha AMS. *Sensor Mater*. 2018;30:179. DOI: 10.18494/SAM.2018.1726.
18. Kolic M, Baglin EK, Titchener S, et al. A 44 channel suprachoroidal retinal prosthesis: Interim functional vision results. *Invest Ophthalmol Vis Sci*. 2020;61:2199.
19. Medical Device Network. Pixium Vision obtains breakthrough device designation for Prima System. 2023. Available at: [www.medicaldevice-network.com/news/pixium-vision-breakthrough-device-prim/?cf-view](http://www.medicaldevice-network.com/news/pixium-vision-breakthrough-device-prim/?cf-view) (Date last accessed: 6 December 2023).
20. Clinicaltrials.gov. Restoration of central vision with the PRIMA system in patients with atrophic AMD (Primavera). Clinicaltrials.gov Identifier: Nct04676854. 2023. Available at: <https://clinicaltrials.gov/study/NCT04676854>.
21. Wang B-Y, Chen ZC, Bhuckory M, et al. Electronic photoreceptors enable prosthetic visual acuity matching the natural resolution in rats. *Nat Commun*. 2022;13:6627. DOI: 10.1038/s41467-022-34353-y.